QA3 Emergency Trolley

Radiotranslucent platform with exceptional X-ray access and low height capability



016/40 05 45 info@dermat.be www.dermat.be



QA3 EMERGENCY TROLLEY SYSTEM: Four new models to choose from






Specially designed for Emergency Departments with radiotranslucent platform, offering exceptional X-ray access as well as a medirail and K8 pressure care mattress.

Emergency Trolley

Description	Code
QA3 Emergency Trolley with Radiolucent Platform	21120



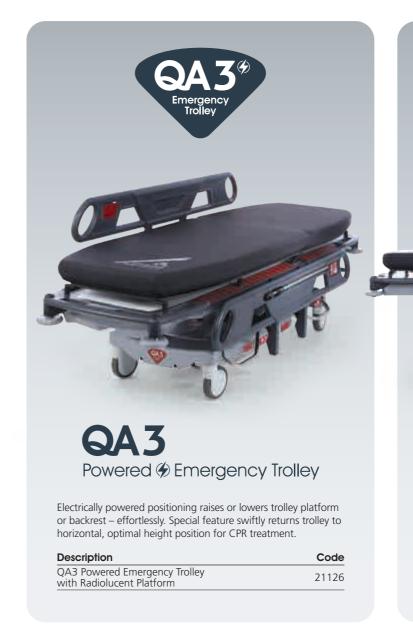



Description	Code
QA3 DRIVE Powered Emergency Trolley with Radiolucent Platform	21124

Quick reference guide to your QA3 Emergency Trolley options: QA3 STANDARD FEATURES Emergency Trolley QA3///IDRIVE STANDARD FEATURES plus INBUILT MOTORISED DRIVE **Emergency Trolley** QA3 **STANDARD FEATURES plus POWERED FUNCTIONS** Powered & Emergency Trolley QA3///IIDRIVE STANDARD FEATURES plus INBUILT MOTORISED DRIVE and POWERED FUNCTIONS Powered & Emergency Trolley

Specialist Bariatric option of any model available: see page 10

Contents

Standard features	4-5
DRIVE features	6-7
Powered features	8-9
Build Options including Bariatric Trolley	10-11
Accessories	12-13
AnetiCare	14-15







Combining industry knowledge with innovative technology



Designed and produced at our UK manufacturing facility using the latest engineering techniques



Manufacturer, designer and supplier of high quality medical equipment around the world



Special bariatric version of this model also available – see page 10







Motorised **IIIIDRIVE** FEATURES











are getting larger, it's crucial that we take moving and handling issues into account, giving practitioners a trolley on which patients can be moved and positioned as easily as possible.

The introduction of our new touch-activated, motorised DRIVE feature offers a major enhancement in a trolley that already has ease of use, practicality and patient comfort built in.

The trolley can also be used in traditional fifth wheel steer mode.

DRIVE handles

innovative touch sensitive design: simply press against the pop-up or fixed yellow pushing handles, and the trolley will move forward; pull back on them, and it will reverse; let go and it will stop



Effortless travel up to gradients of **10**°



Inbuilt motorised drive

no need to push or pull the trolley: our new, inbuilt motorised DRIVE feature linked to the trolley's 5th wheel, offers controlled, effortless travel up to gradients of 10°



Long-life battery – 22 kilometers on a single charge

DRIVE feature powered by advanced nickel-metal hydride rechargeable battery



Highly accurate battery charge indicator –









Motorised DRIVE: effortless travel, ultimate control



Platform height – electrically powered raise and lower mechanism

Offers smooth movement to position the trolley at any point in its height range – 535mm to 835mm (excluding mattress)









Lights on the power panel illuminate to refer the user to the manual if the equipment is being used incorrectly or to indicate when maintenance is required

If a patient needs CPR treatment, there is no time to lose: the powered CPR feature returns the trolley swiftly to a horizontal, optimal height position for **CPR** treatment





Special features for a trolley that's more powerful, by design









Again recognising the trend for increasing patient weights and sizes, we offer special Bariatric versions of the QA3.

The extra wide (800mm) platform significantly increases the surface area for added safety and comfort for Bariatric patients.



capacity

As with all of our trolleys, the QA3 <800 > is fitted with a K8 pressure care mattress.

Specially developed for this Bariatric trolley, the mattress core is made of high density foam, moulded with a unique castellated surface. This increases air flow within the mattress and helps to relieve pressure by distributing the patient's weight as evenly as possible across the surface area – particularly important for heavier patients.

Like other K8 models, the foam has antibacterial and microbial fungistatic properties and its latex-free, 4-way stretch polyurethane (PU) outer cover allows the skin to breathe – for greater patient comfort.

Finally, our unique ultrasonic sealing system effectively welds the seams of the outer cover material together to give an almost unbreakable seal, making the K8 a class leader in infection control.



All models of the QA3 Emergency Trolley have a weight capacity of 320kg and can be fitted with this increased platform width

Description	Code
QA3 Build Option <800> 800mm Extended Width Patient Platform	22160
















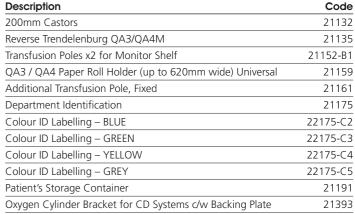











Specialist build options available

AneticAid

(())三田田田

ACCESSORIES

K8 Pressure Care Mattresses

Pressure care and infection control are top priorities in a trolley mattress – and Anetic Aid's independently tested K8 pressure care mattress addresses both of these vital factors.

At its core is high density foam, uniquely specified by Anetic Aid. This means that the mattress can be thinner, while still offering the same, if not better, pressure distribution than a thicker, less sophisticated foam. This is particularly important with tables and trolleys, as it means the overall height can be kept as low as possible for easy access by frail or elderly patients.

A layer of thermally responsive gel-feel memory foam is then added, further relieving pressure by distributing the patient's weight as evenly as possible across the surface area – whatever their position.




Infection control

Both of the foam elements of the mattress have antibacterial and microbial fungistatic properties, while its latex-free, 4-way stretch polyurethane (PU) outer cover is anti-microbacterial, water resistant and vapour permeable. This means it is not only contamination-resistant, but it also allows the skin to breathe – for greater patient comfort.

Finally, a unique ultrasonic sealing system effectively welds the seams of the outer cover material together to give an almost unbreakable seal, making the K8 a class leader in

In the science of Kinetics, there are eight points of contact the human body makes when lying on a hard surface - hence K8.

Description	Code
QA3 K8 Pressure Care Mattress with Zip Cover	21186
QA3 K8 Pressure Care Mattress – Replacement Waterproof Zip Cover	21186-B1

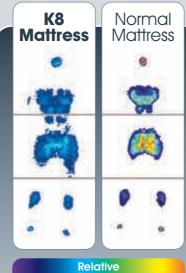

Anetic Aid's K8 technology mattresses can be found in over 90% of NHS hospitals in the UK.

High risk pressure care mattress Independently tested Latex free

Antibacterial/vapour permeable

Welded seams X-ray translucent

Additional products for all QA3 Trolley models

Specially designed and manufactured accessories for QA3 Patient and Emergency Trolley Systems.

Description	Code
Monitor Shelf c/w Removable Refreshment Tray	21152
Refreshment Tray ONLY	21152-A1
QA3 Foot-end Extension c/w 100mm Pad	21153
Theatre Pillow c/w QA3 Mattress Retention Strap	21188
QA3 Foot-end Mattress Protector	21189
QA3 and QA4 Trolley Systems Push / Pull Bar	21390













Emergency Trolley

QA3/IIIIIDRIVE **Emergency Trolley**

Powered **%** Emergency Trolley





Padded Side Rails

Description	Code
QA3 Side Rail Covers, Grey (pair)	21172
QA3 Side Rail Covers, Child Print (pair)	21173





Accessories for radiographers

As well as our Lateral Cassette holder, we offer three additional items to aid patient positioning and help radiographers achieve the best possible image outcomes.





Wallace Support facilitates positioning

of cassette for Modified Axial (Wallace) View patient shoulder X-rays. Consists of two props which attach to Side Rail with spring-loaded guick release handles

Carrier Tray – allows for easier loading and positioning of detectors (cassettes) – while also helping to protect them from accidental damage






Code
21192
21195-A1
21195-B1
21195-C1

Wide range of extra items to complement any QA3 Emergency Trolley







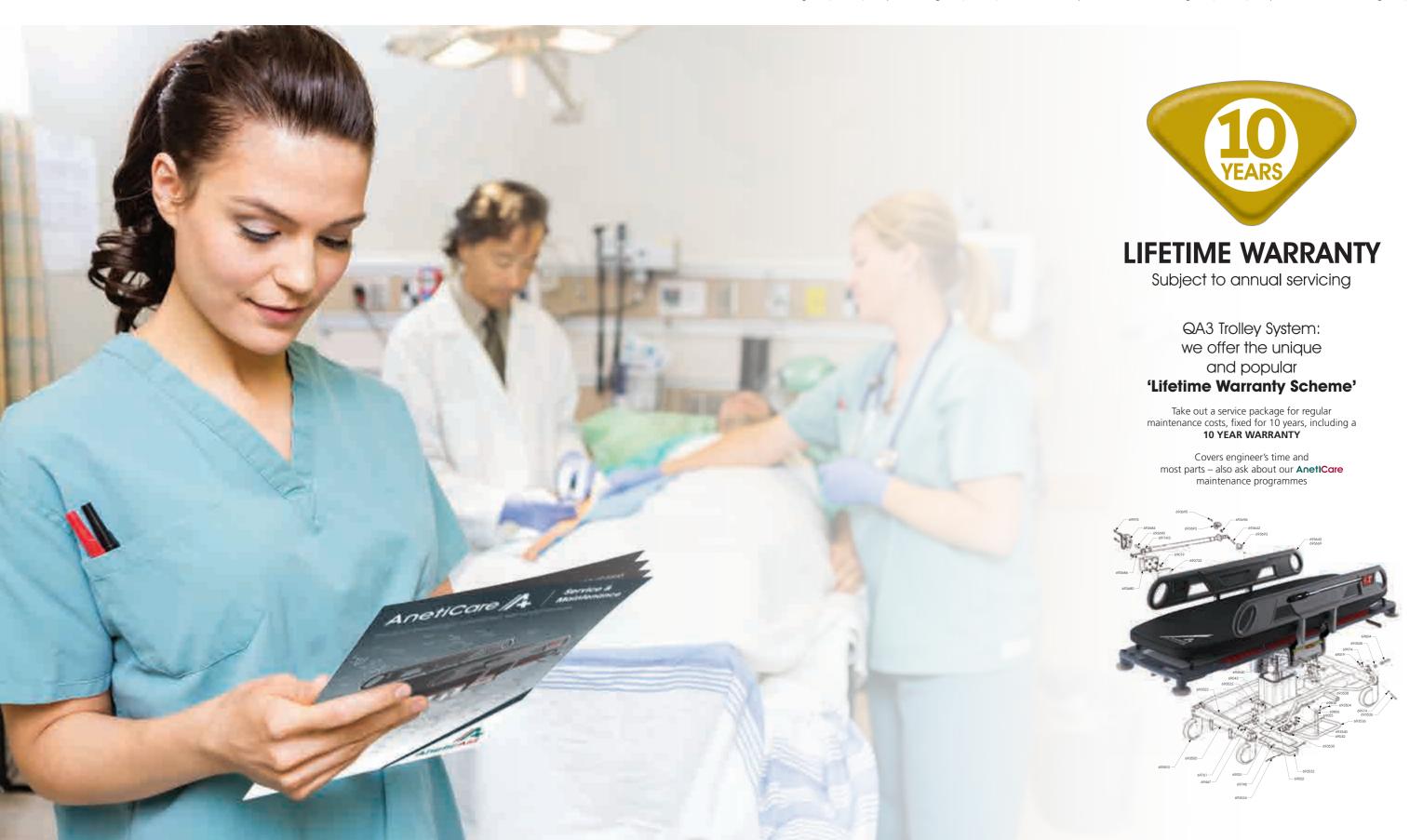








Value for money: minimum downtime, maximum peace of mind



QA3

Emergency Trolley

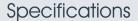

Overall length:

2150mm

Overall width:

835mm (side rails up)

735mm (side rails down)

Height range excluding mattress:

535mm to 835mm

Trolley weight:

148kg (21120)

157kg (21122)

154.5kg (21126)

162.5kg (21124)

Maximum patient weight:

320kg (Please refer to user manual)

Please refer to 'Instructions for Use' for full product specification details

For more information on any of our products or service contracts, call:

T +44 (0)1943 878647

sales@aneticaid.com

aneticaid.com

Innovative medical technology - practically applied

SPECIALISING IN:

Patient & Surgery Trolley Systems

Operation Table Accessories

Stainless Steel Theatre Furniture

Radiation Protection

Electro-Surgical Accessories

Tourniquet Systems

Fibre Optic Instruments

Service & Maintenance

Standard features

X-ray translucent platform – full length radiotranslucent platform, suitable for digital or traditional plates

Hydraulic lift mechanism - exceptional low height capability

Vertically operated side rails - giving a virtually zero transfer gap

Trendelenburg tilt – extremely smooth and controlled mechanism

Latest castor technology – low rolling resistance for high manoeuvrability

5th wheel steering – easily guided by one person

Pushing handles - ergonomically positioned

Gas assisted backlift - horizontal to near vertical

Advanced braking system – exceptionally stable platform

Two fixed transfusion poles – innovative quick release

Cylinder trough – side loading

K8 pressure care mattress – body mould foam, anti-bacterial and latex free

Rotating buffers - at head and foot ends

Infection control – anti-bacterial properties and smooth moulded surfaces









Sapper Jordan Rossi Business Park, Otley Road, Baildon BD17 7AX, UK