



medifa 7000

Mobile
Operating Tables



016/40 05 45 info@dermat.be www.dermat.be

we care.

Advantages at a glance



- With up to three electrical driven joints
- Height adjustment from 580 mm up to 1,100 mm
- Motorised drive by means of fifth castor can be integrated
- Electrical leg rest adjustment, single or parallel
- Optional with longitudinal sliding
- Battery operated back up system can be integrated
- Optimal decubitus prevention due to 80 mm strong upholsteries
- Optimal cleaning and disinfection conditions
- User-friendly control concept
- Suitable for numerous applications and patient positionings
- Suitable for x-ray and c-arm applications
- Wide range of accessories for all surgical disciplines
- Maximum permissible overall load up to 544 kg (1,200 lbs)



Functionality and ergonomics





Motorised drive

- The motorised drive facilitates the positioning of the operating table before and after the surgery (optional, art.-no. 70645)
- Additional handles with a precisely adjustable speed regulator ensure safe and ergonomic patient transport (optional, art.-no. 71570)



Electrical leg rest adjustment

- Enables careful and variable positioning of the patient
- Single or parallel adjustment of leg rests from -90° up to +80°
- Function also applicable for electrical adjustment of back rest in reverse positioning of patient (perfect for shoulder arthroscopy)



Quick and easy handling

 Special coupling points ensure quick and easy locking of the attachment segments and provide an ergonomic workflow



Hygiene and comfort

- Fixing and removing of upholstery elements by means of pin system
- Particularly hygienic and easy to clean (no Velcro fastenings)
- Viscoelastic SF upholstery (H= 80 mm) combines penetration protection and pressure relief
- Electrically conductive and radiolucent

/

Control concept





- 1. Standard: Backlit hand switch with integrated colour display and battery indicator (optional: additional wireless infrared remote control, art.-no. 704070, not retrofittable)
- 2. Standard: Additional control unit at the column
- 3. Optional: Back-up system that ensures functionality of operating table's brake and basic functions in case of an electronic defect (optional, art.-no. 70655, not retrofittable)
- 4. Optional: Foot switch for height, lateral and Trendelenburg adjustment (art.-no. 704060)









 \sim 7

medifa 7000





Mobile Operating Table Art.-no 700200

- Height adjustment from 580 to 1,100 mm (with additional longitudinal sliding from 600 to 1,120 mm)
- Trendelenburg adjustment from -30° up to +30°
- Lateral tilt up to 20°
- Back rest adjustment from -45° up to +80°
- Single or parallel electrical adjustment of leg rests from -90° to +80°
- Possibility of reverse mode
- Flex / Reflex position by hand switch and control unit at column
- Automatic simultaneous leveling of lying surface
- Hand switch and control unit at column for electric adjustments
- Floor locking by 4 extendible stabilization cylinders between the double swivel castors
- Viscoelastic SF upholsteries with memory effect, h = 80 mm, electrically conductive, detachable, radiolucent
- Exchangeable back section segments
- Maximum permissible overall load: 544 kg (1,200 lbs)



Further models:

Art.-no. 700300

Additionally with electrical adjustment of upper back rest and kidney bridge

Art.-no. 71640

Optional: Longitudinal sliding of table top up to 380 mm





- Height adjustment from 580 up to 1,100 mm
- Possibility of patient positioning in reverse mode
- Large X-ray areas in the lying surface guarantee optimal use of C-arms

Configuration example

Art.-no. 700200

Mobile electric operating table

Art.-no. 71640 Longitudinal sliding

Art.-no. 71144 Head rest

Art.-no. 71250_2 Upper back rest

Art.-no. 71615Divided leg rests

Art.-no. 61220SF Arm rest

Price on request!

Specialities

- Bariatric Surgery
- Orthopaedics
- ENT
- Ophthalmology
- Gynaecology
- Urology

Models

Positioning examples







Thoracic and kidney surgery
with head ring, arm rests, body supports and leg
positioning tunnel



with anaesthesia screen, arm restraint cuff, bottom roll and heel pads



Spinal surgery
with rectal positioning device, rectoscopy roll, arm rests
and head positioning pad



Urology and gynaecology with leg rests and head positioning pad

Positioning examples






ENT and Ophthalmology
with foot switch, horseshoe shaped head rest with joint arm,
bottom roll and heel pads

Neurosurgery with bottom roll, arm rests, head positioning pads and

skull clamp system

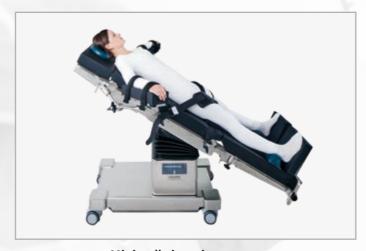

Minimally invasive surgery
with head positioning pad, arm rests, leg straps, heel pads
and foot plates



Gastric banding
with head ring, arm rests and
heel pads



Shoulder arthroscopy
with tripartite back rest, head support, arm rest
and heel pads



Knee arthroscopy
with head ring, arm rest, knee positioning device
and heel pads





Orthopaedicswith extension device and head ring





Arm and hand surgery
with arm extension device and
Weinberger hand cuff


Technical information

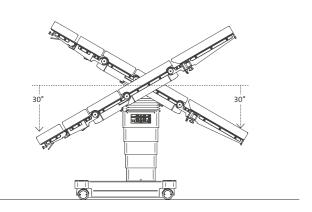


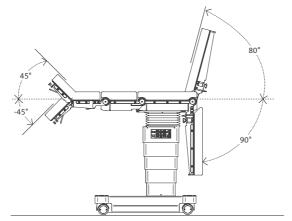


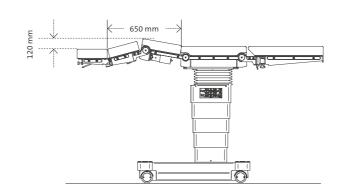

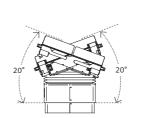


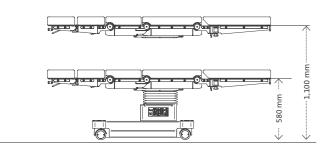

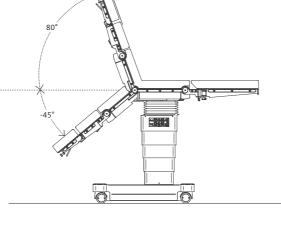


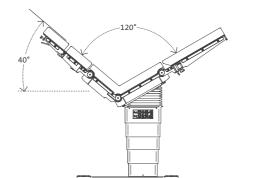




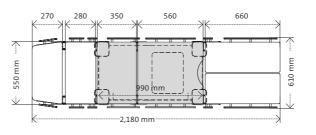


Longitudinal sliding of table top (optional)

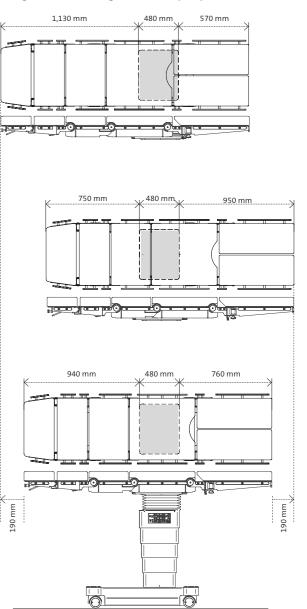

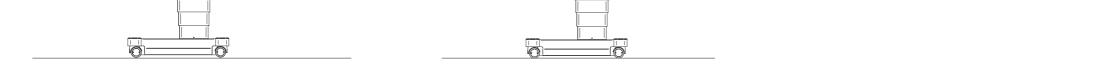



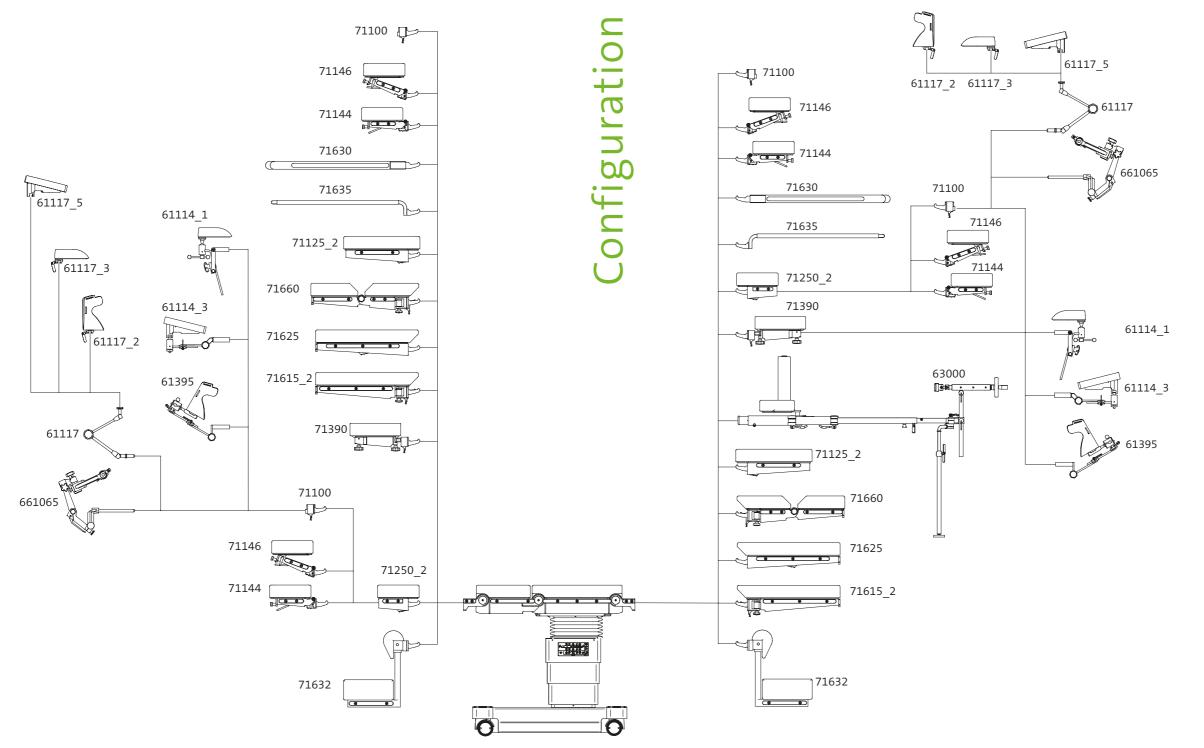



Head rest for medifa 7000 Art.-no. 71144

- Easy adjustment by means of gas spring from -45° to +45°
- Docking and undocking at back rest
- Radiolucent
- L x W: 266 x 550 mm
- Viscoelastic SF upholstery, H = 80 mm, electrically conductive, detachable



Head rest for medifa 7000 Art.-no. 71146

- Easy adjustment by means of gas spring from -45° to +45°
- Double joint enables additional tilting
- Docking and undocking at back rest
- Radiolucent
- L x W: 266 x 550 mm
- Viscoelastic SF upholstery, H = 80 mm, electrically conductive, detachable



Head-end rail Art.-no. 71148

- Additional rail at head rest
- Length: 300 mm
- Material: Stainless steel
- Suitable for: 71144, 71146
- Piece



Side rails at head rest Art.-no. 71149

- Additional side rails at both sides of head rest
- · Length: 160 mm
- Material: Stainless steel
- Suitable for: 71144, 71146
- Pair



Upper back rest for medifa 7000 Art.-no. 71250 2

- Docking and undocking at lower back rest
- Radiolucent
- L x W: 268 x 550 mm
- Viscoelastic SF upholstery, H = 80 mm, electrically conductive, detachable
- Integrated side rails (EU)



Tripartite back rest for medifa 7000 Art.-no. 71390

- Incl. head rest
- Easy access to shoulder area by means of three detachable segments
- Docking and undocking at lower back rest
- Radiolucent
- L x W: 320 x 550 mm
- Viscoelastic SF upholsteries, H = 80 mm, electrically conductive, detachable



Pelvis section extension for medifa 7000 Art.-no. 71125 2

- Extension of pelvis section by 400 mm
- With cut-out
- Radiolucent
- L x W: 400 x 550 mm
- Viscoelastic SF upholstery, H = 80 mm, electrically conductive, detachable
- Integrated side rails (EU)



Leg rest for medifa 7000

Art.-no. 71625

- Electrically adjustable
- Docking and undocking at pelvis section
- Radiolucent
- L x W: 660 x 550 mm
- Viscoelastic SF upholstery, H = 80 mm, electrically conductive, detachable
- Integrated side rails (EU)



Divided leg rest for medifa 7000 Art.-no. 71615_2

- Electrically adjustable
- Swivelling up to 180°
- Docking and undocking at pelvis section
- Radiolucent
- L x W: each 660 x 266 mm
- Viscoelastic SF upholsteries, H = 80 mm, electrically conductive, detachable
- Integrated side rails (EU)
- Pair





Double divided leg rest for medifa 7000 Art.-no. 71660

- Electrically adjustable
- Swivelling up to 180° and bendable
- Docking and undocking at pelvis section
- Radiolucent
- L x W: each 662 x 266 mm
- Viscoelastic SF upholsteries, H = 80 mm, electrically conductive, detachable
- Integrated side rails (EU)
- Pair

Carbon positioning rest for medifa 7000 Art.-no. 71630

- Incl. upholstery
- For intraoperative x-ray, optimal radiolucency (360°)
- Docking and undocking at pelvis section or lower back rest
- L x W: 750 x 540 mm
- Viscoelastic SF upholstery, H = 70 mm, electrically conductive, detachable

Transfer leg rest for medifa 7000 Art.-no. 71635

- For short-time positioning of legs
- Docking and undocking at pelvis section
- Radiolucent
- L x W: 900 x 500 mm





Extension segments Art.-no. 61355SF

- For extension of the operating table's lying surface
- Integrated side rails (EU)
- L x W: each 310 x 110 mm
- Viscoelastic SF upholsteries, H = 70 mm, electrically conductive, detachable
- Suitable for all operating tables with viscoelastic foam upholsteries
- Pair



Foot switch for medifa 7000 Art.-no. 704060

- Additional control unit to adjust the main functions height, lateral and Trendelenburg movement
- Tethered



Infrared remote control for medifa 7000 Art.-no. 704070

- Wireless data transfer
- Control panel for all electric functions of operating table
- Incl. charging station



Longitudinal sliding for medifa 7000 Art.-no. 71640

- Electrical longitudinal sliding of table top up to 380 mm
- Changes lowest height of operating table from 580 mm to 600 mm
- Cannot be retrofitted!



Motorised drive Art.-no. 70645

- Motorised fifth castor for manoeuvring and positioning of operating table
- Cannot be retrofitted!



Hand grips for motorised drive Art.-no. 71570

- For manoeuvring and positioning of operating table by means of electronic drive at operating table's head end
- With integrated speed control
- Material: Stainless steel
- Suitable for: 70645



Back-up system for medifa 7000 Art.-no. 70655

- System ensures the operability of operating table's basic functions (height, lateral, Trendelenburg and brakes) in case of electronic defect by means of 24V battery
- Cannot be retrofitted!



Radiographic cassette carrier for medifa 7000 Art.-no. 71643

- To insert radiographic cassettes under table top of operating table
- With telescopic rod
- L x W x H: 470 x 385 x 16 mm
- Suitable for radiographic cassettes with dimensions 14" x 17" (US)
- Material: Stainless steel
- Suitable for: 71642_2, 71642_3





Radiographic cassette tunnel for 700200 Art.-no. 71642_2

• To insert radiographic cassette carrier

Radiographic cassette tunnel for 700300 Art.-no. 71642_3

• To insert radiographic cassette carrier



Additional accessories

A wide range of optional accessories complements the application spectrum of our operating tables.

For central storage of accessories we offer a trolley (art.-no. 341045) made of polished stainless steel.

Further accessories can be found in our brochure "Accessories for operating tables".





medifa



The medifa healthcare group unites competence concerning operating rooms under one roof.

medifa is an owner-run and internationally active complete solution provider in the medical technology.

The product and service portfolio covers all operating room requirements from modular wall, door, ceiling and cabinet systems (RooSy® Room Systems) for the turnkey operating room, as well as mobile operating tables, operating table accessories and operating lights, right up to examination chairs and medical

furniture. In addition, one of the largest German OEM suppliers in the medical technology sector also belongs to medifa.

medifa GmbH & Co. KG (Finnentrop), medifa hygienic rooms GmbH (Ötigheim) and medifa metall und medizintechnik GmbH in Rastatt are all part of this group of companies. The medifa healthcare group employs approximately 350 employees at three locations within Germany. The planning, engineering and production


take place exclusively in Germany, and with over 80% added value in the group of companies, the group guarantees fast reaction times.

The slogan "we care" means for us that we consistently orientate ourselves according to the requirements of our customers. Customers all over the world can rely on top-class consultation, processing and quality from a single source, made in Germany according to certified national and international quality standards.

Products

- Modular operating room systems (RooSy)
- Mobile operating tables and accessories
- Operating lights (H!Light)
- Examination and treatment chairs
- Medical furniture and trolleys
- Customized medical solutions

medifa worldwide

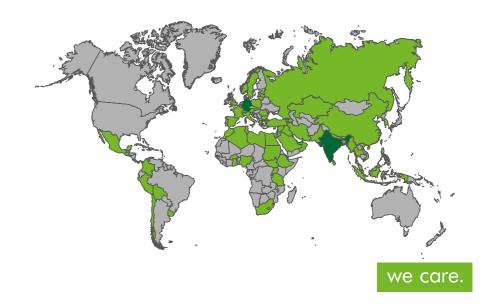













Version 10-20-EN

Subject to modifications due to technical development. This document is intended to provide information to an international audience outside of the US.



medifa GmbH & Co. KG

Industriestraße 5 57413 Finnentrop Germany fon +49 2721 71 77 -0 fax +49 2721 71 77 -255 info@medifa.com

we care.