

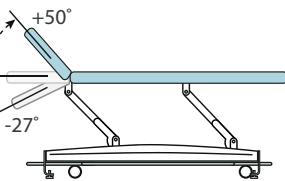
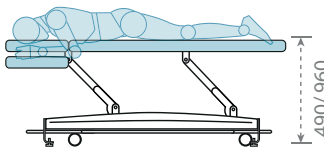
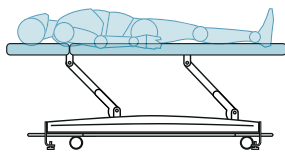
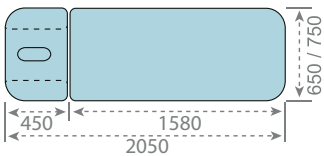
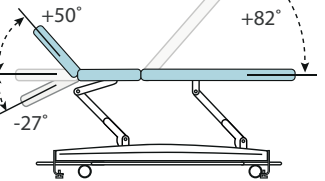
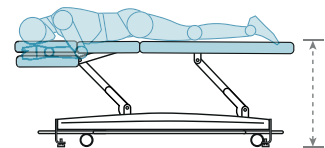


Technical data





Model: **MTS2-H/E-65/75**

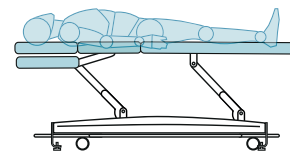





Model: **MTS3-H/E-65/75**

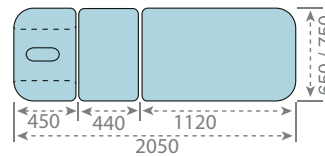





Model: **TTS2-H/E-65/75**

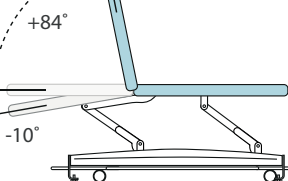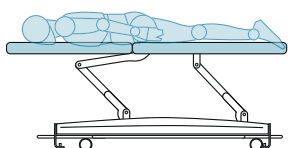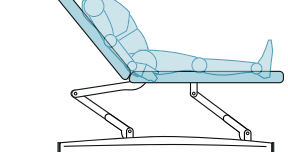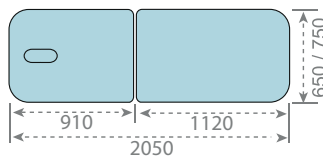





Model: **TTS3-H/E-65/75**

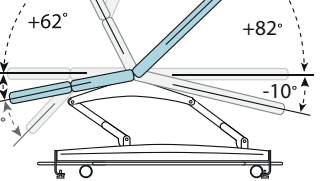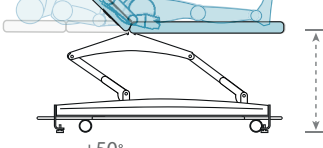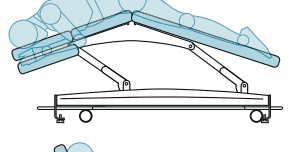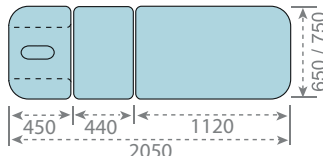



Model: **TTS4-65/75**

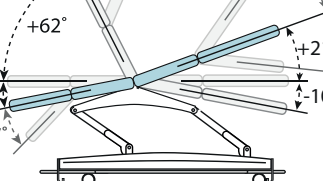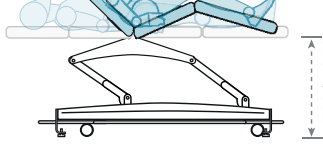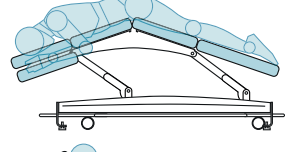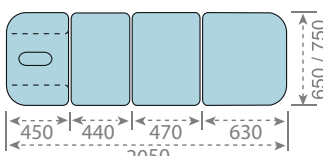





Model: **TTS8-65/75**

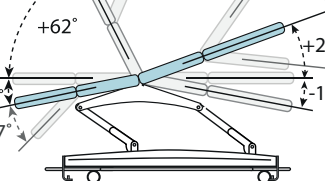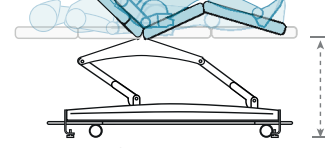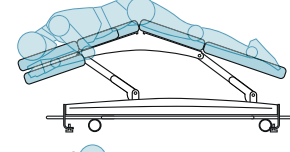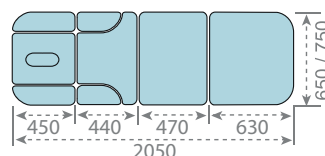

Model: **BTS-100/120**

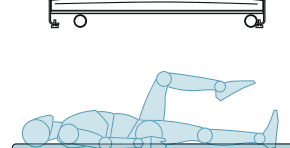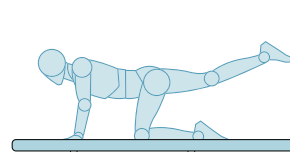

Model: **BTS2-100/120**

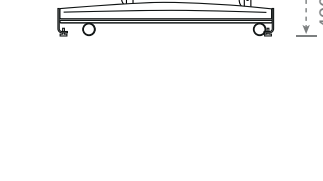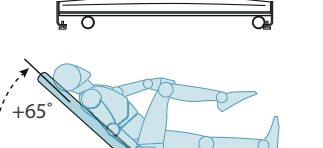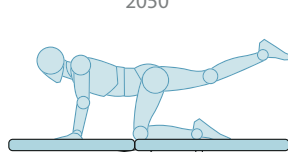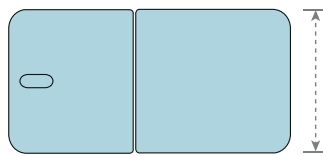

Equipment and accessories

standard ● optional ○

	MTS2-E	MTS2-H	MTS3-E	MTS3-H	TTS2-E	TTS2-H	TTS3-E	TTS3-H	TTS4	TTS8	BTS	BTS2
durable, soft foam upholstery (8 cm) – SER-SFU8-S	●	○	●	○	●	○	●	○	●	●	●	●
durable, firm foam upholstery (8 cm) – FFU8-S	○	○	○	○	○	○	○	○	○	○	○	○
breathing hole with welded and taped seams (plug included)	●	●	●	●	●	●	●	●	●	●	●	●
synchronously adjustable three-part head section – TVL-S	○	○	○	○	/	○	○	○	○	●	/	/
paper roll holder under head section (not available with TVL-S) – NZP1-S	○	○	○	○	○	○	○	○	○	/	/	/
paper roll holder at table base – NZP2-S	○	○	○	○	○	○	○	○	○	○	/	/
PVC cover for leg section – PVC8-S	○	○	○	○	○	○	○	○	○	○	○	○
retractable castors with soft release (Ø 75 mm) – SER-RCS	●	●	●	●	●	●	●	●	●	●	●	●
hand switch	/	/	/	/	/	/	/	/	●	●	/	/
circular switch	/	/	●	/	●	/	/	/	●	●	/	/
safety switch – EVV-S	○	/	○	/	○	/	○	/	○	○	○	○
battery backup – BB-S	○	/	○	/	○	/	○	/	○	○	○	○
plug connector for potential equalization – PPEQ	○	○	○	○	○	○	○	○	○	○	○	○
safe working load 250 kg	●	●	●	●	●	●	●	●	●	●	●	●



reddot winner 2020



Dermat⁺

info@dermat.be
www.dermat.be

MASSAGE AND THERAPY TABLES S

We make it our priority that the treatment table fulfills the specific requirements of your practice. Visit our **ON-LINE CONFIGURATOR** and equip the table so that it perfectly supports your daily tasks and provides your patients with safety and maximum comfort.



CONFIGURE



Head section

- breathing hole with memory foam and matching plug
- welded and taped seams around the breathing hole preventing the penetration of massage oils and lotions
- practitioner's unhindered access under three-part head section
- a one-handed push down or up adjusts both armrests at once

Upholstery

- durable and seamless
- soft or firm double-layer foam padding, 8 cm thick
- wider lying surface due to the padding finished at a right angle

Hand switch includes CPR button for fully extended position, SOS button for Trendelenburg position and a memory function to store a position of practitioner's choice.



Table base

- ABS plastic cover is screwless for easier cleaning
- standard 250 kg safe working load

Stronger control unit

- new software enables synchronous movement of table's middle two sections
- easier to replace

Soft retractable system

- includes a shock absorber to soften the release of castors

Circular switch

- simple up & down operation from all sides of the table

So the cleaning becomes easier

To ensure **simple cleaning** and **high hygiene standards**, the **UPHOLSTERY** is **SEAMLESS** and made of **high-quality BIOCOMPATIBLE imitation leather**. **SEAMS** around the breathing hole are **WELDED** and **TAPED** to **prevent penetration** of any **liquids, massage lotions and oils** in particular. The **ABS plastic COVER** that protects the table base is **screwless** and thus completely **smooth**. The head and leg section **ADJUSTMENT LEVERS** with a mechanism **integrated** in the frame are **easier to press** and **clean**.

To maximize the comfort and ease the use

To maximize the patient comfort, the **8 cm thick PADDING** is made of **two-layer foam**. It is **durable and comfortable** along the **whole lying surface** which is now even wider due to the **padding finished at a right angle**. To fully support the treatments, you can choose between **soft** and **firm** foam padding, which is available in more than 30 colors. The **BREATHING HOLE** with **memory foam** on the edges ensures perfect comfort in prone position. A **matching PLUG** is standard too. Optionally, tables can include a **THREE-PART HEAD SECTION** that allows **practitioners** an **unhindered access** when they treat patients seated at the head section. With the **NEW MECHANISM** for armrest adjustment, only a **push by one hand** is enough to lower or raise both armrests at once.

You can choose between **ELECTRICALLY** or **HYDRAULICALLY ADJUSTABLE** examination tables. Simple height adjustment enables **easy access** for patients as well as **comfortable working height** for the practitioner. By pushing the **CIRCULAR SWITCH up** or **down** you can **lower** or **raise** electric tables completely effortlessly from **all four sides** of the table. The height of the hydraulic tables can be adjusted by pressing the **FOOT PEDAL** on **two sides** of the table. The **HAND SWITCH**, beside **adjusting the height and tables' sections**, also includes two buttons for the **prestored completely flat and Trendelenburg positions** and a **memory function** to store a position of practitioner's choice.

The standard **SAFE WORKING LOAD** has been increased to **250 kg**, which contributes to the **perfect stability**. Practitioners can rely on the table even if they apply strong pressure during treatments. Due to the **NEW SOFTWARE** that enables **synchronous movement** of the middle two lying sections, they can **position** their **patients easier** and **quicker**. The table stands firmly on its **RUBBER FEET**. To move the table across the room, it features **RETRACTABLE CASTORS**. Due to the **shock absorber**, the release of the castors is softened.

For various treatments, massages and therapies

The **BREATHING HOLE** with a **MATCHING PLUG** and **ELECTRIC HEIGHT ADJUSTMENT** are now the **standard equipment** of the renewed S-line massage and therapy tables. Two massage table models offer a **wide and comfortable lying surface** to perform massages and therapies. If treatments demand positions with negative section incline to allow the **Jackknife** or **Trendelenburg position**, 3- , 4- and 8-section therapy table models are available. The **fully equipped** Therapy table S8 standardly includes a **three-part head section** and **sturdier foldable side wings** integrated in the back section. Bobath table offers **WIDE SURFACE** for rehabilitation treatments. Now both, one and two-section model, include a **HAND SWITCH** to adjust the height and with a two-section Bobath the head section as well.